

xylem



# OPTi<sup>®</sup> Digital Handheld Refractometer

THE PORTABLE REFRACTOMETER FOR MEASUREMENTS ON-THE-GO



a xylem brand

Portable. Reliable.  
Digital. OPTi®.



# Portable. Reliable. Digital. OPTi®.

The pocket sized digital refractometer that's ready when you are.

The OPTi Digital Handheld Refractometer from Bellingham + Stanley (pioneers of the handheld refractometer) is a highly reliable, premium quality instrument used to control the dissolved solids or blend ratios of products ranging from fresh fruit to industrial chemicals.

## Portable and precise

Portable and simple to operate, the OPTi provides excellent "on-the-spot" concentration analysis with automatic temperature correction. Accurate one button measurement and a clear display presents reliable results, no matter who takes the reading.

Thanks to its compact design and clever marriage of smart software and ergonomic hardware, OPTi refractometers are perfect for taking measurements anywhere.



## Software

- Quick 2 second read time of sample
- Onboard library of over 40 scales
- Up to 3 active scales on a single device
- High Ambient Light (HAL) detection warns of too much light entering the sensor
- Unique "AG Test Mode" that allows the use of long-life non sucrose based certified reference materials
- Programmable read delay for sample temperature stabilisation

## Hardware

- Clear digital LCD display
- Ergonomic & durable handheld design
- IP65 rated
- Easy-clean, stainless steel prism dish
- Rapid temperature stabilisation

A high performance refractometer with a vast choice of scales for a wide range of applications – all in the palm of your hand.

## Bellingham + Stanley gives you more from your refractometer

The OPTi revolutionises the way digital handheld refractometers are supplied and purchased. With the OPTi refractometer you will have all of the scales you'll need, right in the palm of your hand – no paying extra for high brix models.

An onboard library of over 40 common scales including Brix, Refractive Index, °Baumé, °Butyro, Colostrum Quality, Wort, Ethylene Glycol % and many more can be swapped and selected at any time.

### As easy as 1, 2, 3

The choice of scale is yours. With up to 3 scales readily available at the touch of a button, the OPTi is the most versatile and feature rich refractometer on the market.

Need Brix for one application, and Alcohol Probable for the next? No problem! Switching scales only takes a moment, and with a read time of 2 seconds per sample, the OPTi is one of the quickest and most convenient ways to monitor and control your concentration on-the-go.

### Full 0-95 Brix Scale

The Brix scale is the most common scale used in refractometry, yet many refractometers cover only a small portion of it. OPTi gives you more, with coverage of the full Brix range (0-95) at no extra cost.

### Over 40 Scales

Each OPTi comes with 3 interchangeable scales pre-loaded from an onboard library of over 40 individual choices - meaning there are over 9880 possible OPTi combinations.

### Versatile


Thanks to its durable design and vast range of scales, the OPTi is suitable for a wide range of applications; from the clinical laboratory to the harshest of environments.

### Digital Display

Readings from an OPTi refractometer are clear and defined thanks to the bold digital LCD display.

### Ready when you are

The refractometer that fits in your pocket. With a long battery life, true portability, and durable design: OPTi refractometers are ready when you are.

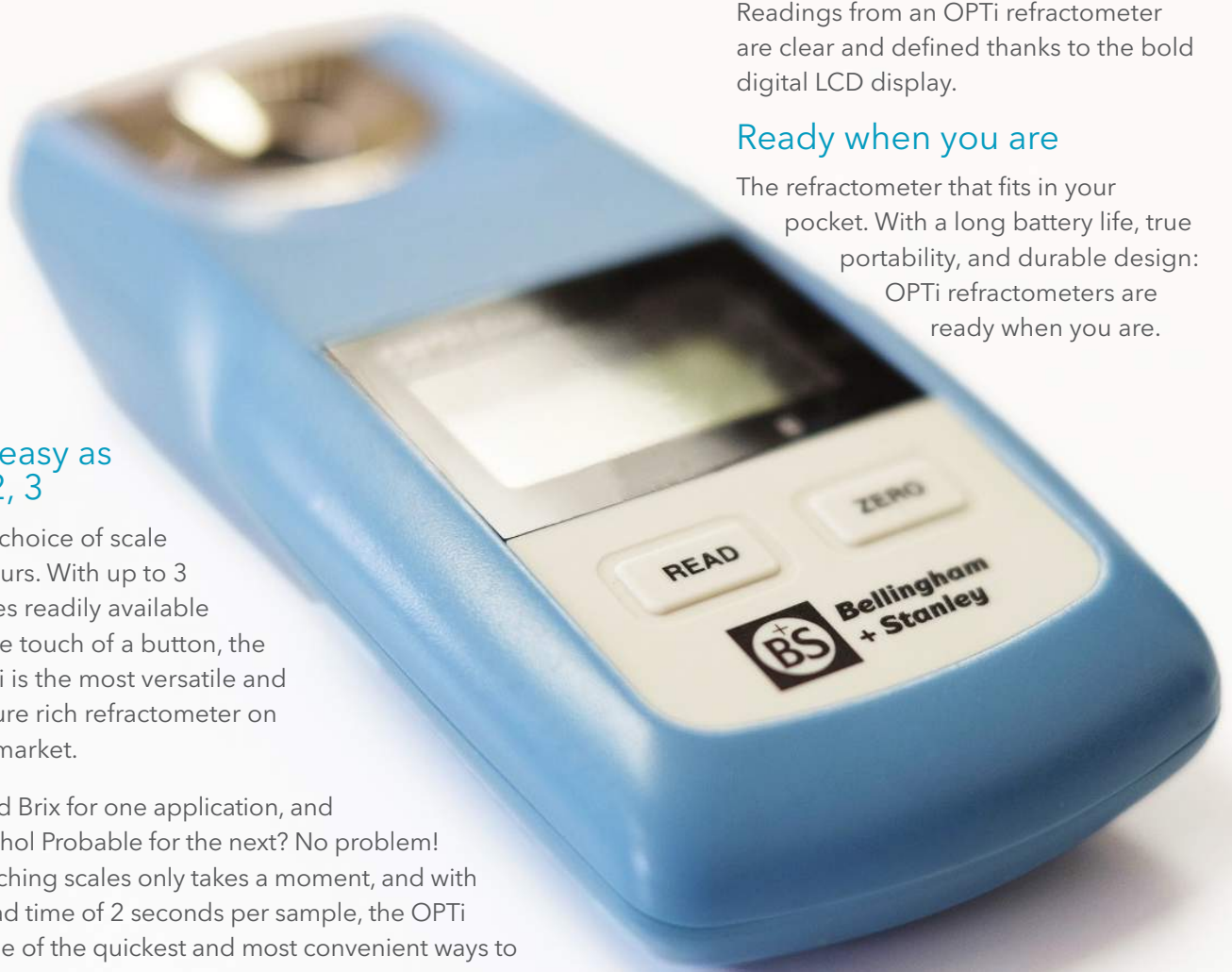

# An OPTi® For All Occasions.

## The convenient way to measure and control on-the-go

The OPTi Digital Handheld Refractometer can be used for measuring fruit ripeness, fruit juices and concentrates, carbonated beverages, alcohol content and grape must, confectionery and jams, sugar and many other food applications. Industrial applications include the measurement of soluble oils, quenchants, glycols, antifreeze, air-conditioning and heat transfer fluids, aviation fuel inhibitors, aeroplane wing de-icing surfactants, starch, textile surfactants and more.



### Food & Beverage

Used for testing the ripeness of fresh fruits such as grapes and tomatoes simply by squeezing a small amount of juice directly onto the refractometer prism. OPTi refractometers are used in the manufacturing of jam, marmalade, syrup and other high sugar content products as well as quickly testing for moisture in honey direct at the hive. They are ideal for testing finished juices and other soft drinks. **Scales include:** Brix, °Butyro, HFCS, Salinity, Water in Honey, and more.



### Beer & Wine

From home-brew beer to the microbrewery and beyond, our refractometers have been helping brewers for years. In the beer world, refractometers are used to measure Wort prior to fermentation as well as determining alcohol content. We supply OPTi refractometers to vineyards around the globe so that grape ripeness can be checked before harvest, sugar content assessed and alcohol content measured in the final stages of production. **Scales include:** Brix, Wort SG, °Baume, Alcohol Probable, Oechsle, and more.



### Industrial

Bellingham + Stanley's industrial scales cover a broad spectrum of applications. A common scale used within the industrial sector is Glycol % - useful for engineers working on air-conditioners, heat exchangers or even pasteurisers. When using the Brix scale, refractometers can be used for fire fighting foam and industrial coolants. In the aviation industry there is a need for determining the concentration of aircraft anti-icing fluids.

**Scales include:** Brix, RI, Calcium Chloride %, Ethanol %, Ethylene Glycol %, and more.



## Automotive

OPTi refractometers are able to test coolant fluids (measurable in either °F or °C frost protection) and battery acid, as well as AdBlue® (DEF) concentration - ideal for vehicle service centres, garages, car manufacturers and fleet managers. Keep your vehicles running smoothly and safely with OPTi. **Scales include:** Brix, RI, Sulphuric Acid SG, AdBlue®, and more.



## Veterinary & Life Science

Bellingham + Stanley offer a selection of scales for veterinarians and farmers for testing colostrum, blood and urine samples taken from animals. Thanks to needing only a small drop of sample to give a clear and accurate reading, OPTi refractometers are perfect for measuring samples that are often difficult to obtain or, in the case of colostrum, being a valuable commodity. **Scales include:** Brix, RI, Colostrum Quality, Urine Specific Gravity, Serum Protein, and more.



## And Many More...

There are many more applications where OPTi digital refractometers can help you measure and control, and as such Bellingham + Stanley has developed over 40 different scales to support your business. Visit our website or speak to the Bellingham + Stanley Customer Care team today.

## All the accuracy you need in an easy-to-use handheld device

All OPTi refractometers feature a four decimal place LCD display, Temperature Compensation (ATC), robust design with IP65 rating, a choice of three scales in a single device.

Typical Specifications	Sucrose (Brix)	Refractive Index	Specific Gravity (SG)	Concentration %
<b>Range</b>	0-95	1.33-1.53	1.000-1.050	0-95
<b>Resolution</b>	0.1	0.0001	0.001	0.1
<b>Accuracy</b>	±0.2	±0.0003	±0.003	±0.2
<b>Temperature Compensation (ATC)</b>	Brix (ICUMSA) Application Specific or Brix			
<b>Temperature Range</b>	5-40°C (Ambient)		5-80°C (Sample)	

For a full list of available scales please visit [www.bellinghamandstanley.com](http://www.bellinghamandstanley.com)



# Precision analysis in your pocket.

## True Portability

Fits in the palm of your hand and will give you over 10,000 readings on 2 x AAA batteries.

## Stainless Steel Prism Dish

Rapidly stabilises sample temperature.

## Durable Design

Made from ultrasonically welded ABS plastic, making it resistant to physical bumps, chemical corrosion, and water ingress to IP65.

## Clear Digital Display

The clear LCD display helps you record readings quickly and easily.

## Calibration

Requires only water to perform a ZERO calibration.

## Simple Operation

Apply the sample to the prism dish and press READ.



# Xylem |'zīləm|

- 1) The tissue in plants that brings water upward from the roots;
- 2) a leading global water technology company.

Bellingham + Stanley is part of Xylem Analytics and a leading provider of refractometers and polarimeters.

Xylem Analytics' global brands have been leaders in the laboratory instrumentation market for decades, and are relied upon every day across more than 150 countries. Working in true partnership with our clients, we listen, learn and adapt to individual needs, offering deep application expertise built upon our long history of innovation in instruments and services. Our solutions for analysis, measurement and monitoring help enable many of today's modern laboratories and industrial processes, and provide our customers the trusted and high performing solutions they need to succeed.

Xylem Analytics is part of Xylem Inc., a global company focused on solving the world's most challenging and fundamental water issues. As accurate analysis is crucial to the water industry, Xylem Analytics taps its diverse product brands for leadership in that field and beyond, providing the best laboratory and field monitoring instrumentation across a wide variety of industries.

**For more information on how Xylem can help you, go to [www.xylem.com](http://www.xylem.com)**



**Unit 2, 11 Narloo Street, Malaga, Western Australia 6090**  
**Tel: +618 6240 5600 Fax: +618 9209 3944**  
**E-mail: [sales@perthscientific.com.au](mailto:sales@perthscientific.com.au)**